

Directions for Use

PD Tilt Urine Drug Test Cup (AS/NZS 4308:2023)



Prepare testing consumables.
Check expiry date before use.




Provide cup to donor to provide a urine sample above the minimum urine level.




After urine sample is provided, screw lid on tightly.




Check temperature by reading the green colour indicator (must be between 32°C and 38°C).

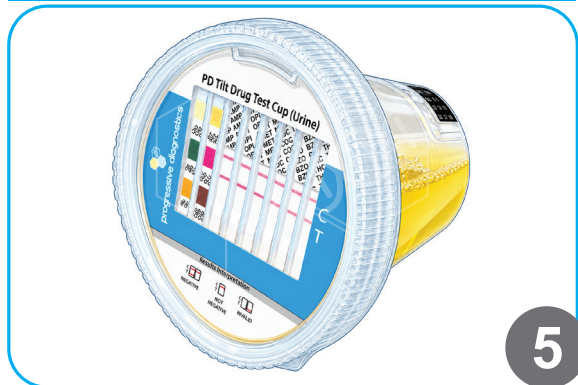


Tilt cup on its side to activate the test.

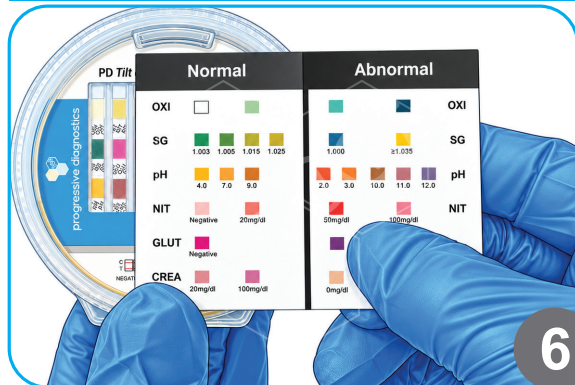

Read adulteration tests. First strip (OXI, SG, PH) and Second strip (NIT, GLUT, CREA).

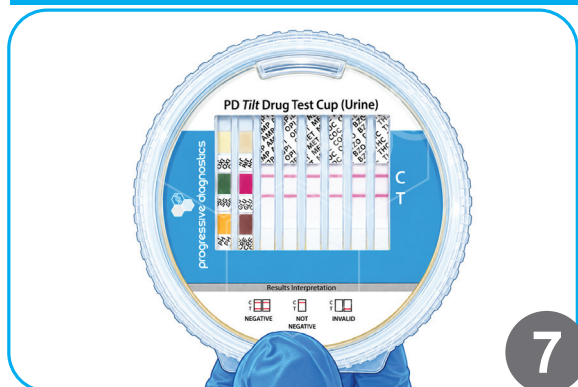

Read drug test strips. If test line(s) not appearing wait at least 5 minutes before result interpretation.

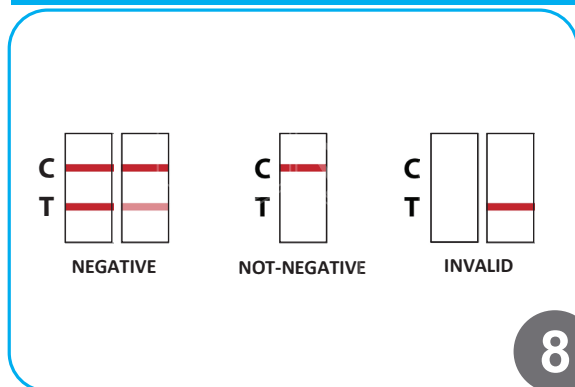

Result interpretation.

